



1 Methodology

Antibiofilm activity was quantified using crystal violet-stained biofilm biomass (OD570) only. For each strain/condition, blank-corrected OD570 was compared against the blank-corrected untreated Positive control (matched by media condition): % Inhibition = $100 - [(OD570 \text{ sample} - \text{blank}) / (OD570 \text{ positive control} - \text{blank}) \times 100]$. No OD600 growth normalization was used for this metric.

Classification thresholds applied — Weak inhibition: $\geq 25\%$; Moderate inhibition: $\geq 50\%$; Strong inhibition: $\geq 75\%$; below the weak threshold, wells were classed No/low inhibition.

2 Instrument & acquisition details

OD570 file — Device: — Wavelength: nm Date: 08 May 2026 Time: 12:15:51 PM Temp: —

3 Summary statistics

Strains/Conditions analysed: 20 Media conditions: 1 Blank wells: 3

Strong inhibition: 17 Moderate inhibition: 1 Weak inhibition: 0 No/low inhibition: 2

4 % Inhibition by strain/condition — Figure 1



BiofilmQ-HT · Antibiofilm Activity Report

% Biofilm Inhibition · crystal violet (OD570) vs. untreated positive control

Generated 08 Aug 2026

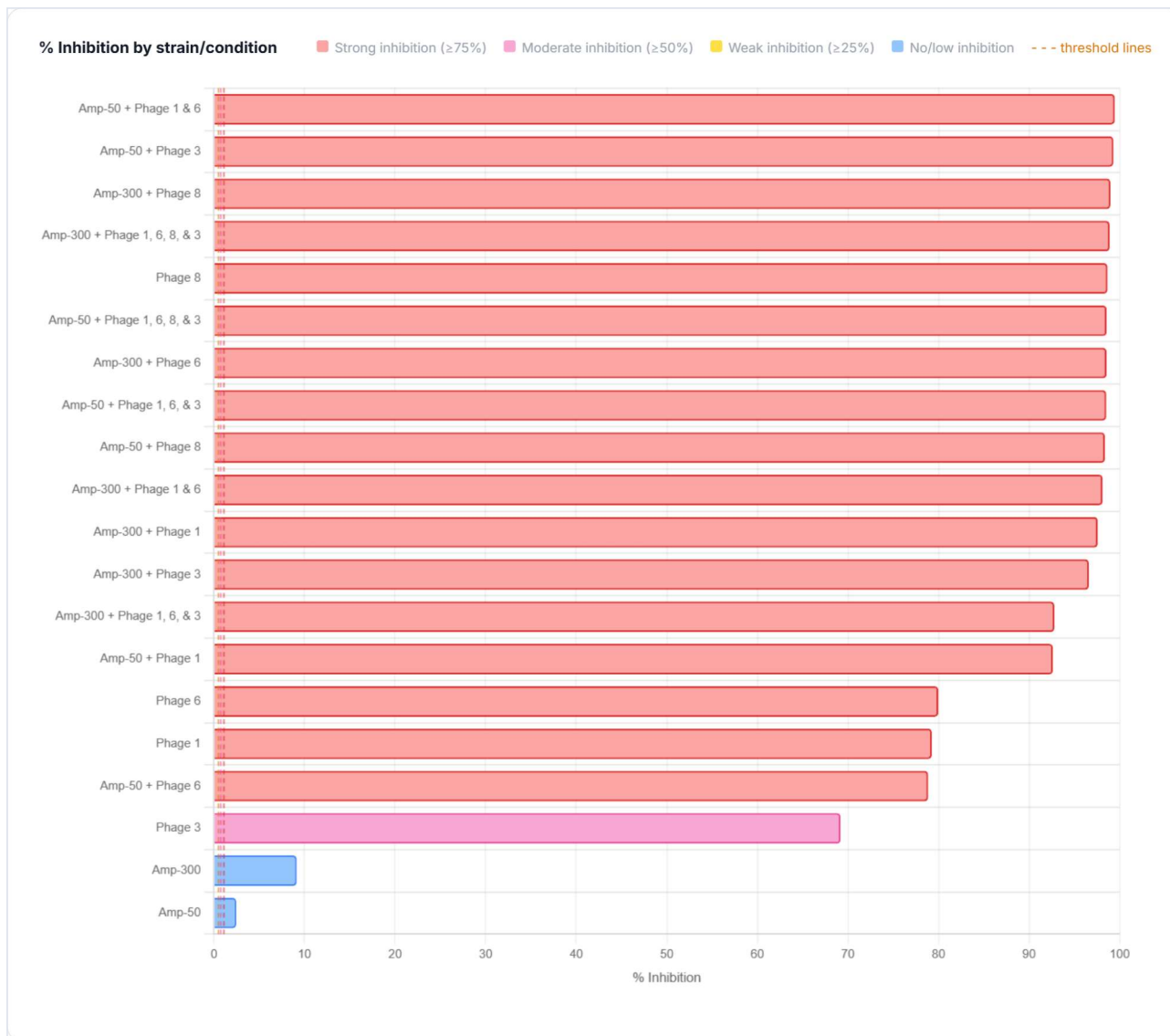




Figure 1. % Inhibition per strain/condition, colour-coded by classification. Dashed lines mark the weak/moderate/strong thresholds.

5 Detailed results table — Table 1

Strain/Condition / ID	Type	Media	n	Mean OD570	Blank OD570	% Inhibition	Classification
Amp-50 + Phage 1 & 6	Strain/Condition	Default	3	0.122	0.115	99.5%	Strong inhibition
Amp-50 + Phage 3	Strain/Condition	Default	3	0.124	0.115	99.3%	Strong inhibition
Amp-300 + Phage 8	Strain/Condition	Default	3	0.128	0.115	99.0%	Strong inhibition
Amp-300 + Phage 1, 6, 8, & 3	Strain/Condition	Default	3	0.129	0.115	98.9%	Strong inhibition
Phage 8	Strain/Condition	Default	3	0.132	0.115	98.7%	Strong inhibition
Amp-50 + Phage 1, 6, 8, & 3	Strain/Condition	Default	3	0.133	0.115	98.6%	Strong inhibition
Amp-300 + Phage 6	Strain/Condition	Default	3	0.133	0.115	98.6%	Strong inhibition
Amp-50 + Phage 1, 6, & 3	Strain/Condition	Default	3	0.133	0.115	98.5%	Strong inhibition
Amp-50 + Phage 8	Strain/Condition	Default	3	0.135	0.115	98.4%	Strong inhibition
Amp-300 + Phage 1 & 6	Strain/Condition	Default	3	0.138	0.115	98.1%	Strong inhibition
Amp-300 + Phage 1	Strain/Condition	Default	3	0.145	0.115	97.6%	Strong inhibition
Amp-300 + Phage 3	Strain/Condition	Default	3	0.157	0.115	96.6%	Strong inhibition
Amp-300 + Phage 1, 6, & 3	Strain/Condition	Default	3	0.205	0.115	92.8%	Strong inhibition
Amp-50 + Phage 1	Strain/Condition	Default	3	0.207	0.115	92.6%	Strong inhibition
Phage 6	Strain/Condition	Default	3	0.365	0.115	80.0%	Strong inhibition
Phage 1	Strain/Condition	Default	3	0.374	0.115	79.3%	Strong inhibition
Amp-50 + Phage 6	Strain/Condition	Default	3	0.378	0.115	78.9%	Strong inhibition
Phage 3	Strain/Condition	Default	3	0.499	0.115	69.2%	Moderate inhibition
Amp-300	Strain/Condition	Default	3	1.248	0.115	9.1%	No/low inhibition



Strain/Condition / ID	Type	Media	n	Mean OD570	Blank OD570	% Inhibition	Classification
Amp-50	Strain/Condition	Default	3	1.331	0.115	2.5%	No/low inhibition
K. pneumoniae	Positive ctrl	Default	3	1.362	0.115	—	Pos ctrl

Table 1. Per-strain/condition results, blank- and media-matched. n = number of replicate wells averaged.

6 Plate heatmap — Figure 2





BiofilmQ-HT · Antibiofilm Activity Report

% Biofilm Inhibition · crystal violet (OD570) vs. untreated positive control

Generated 08 Aug 2026

Figure 2. 96-well plate layout coloured by classification, with per-well % inhibition shown.